




Proudly made  
and manufactured  
in Canada!

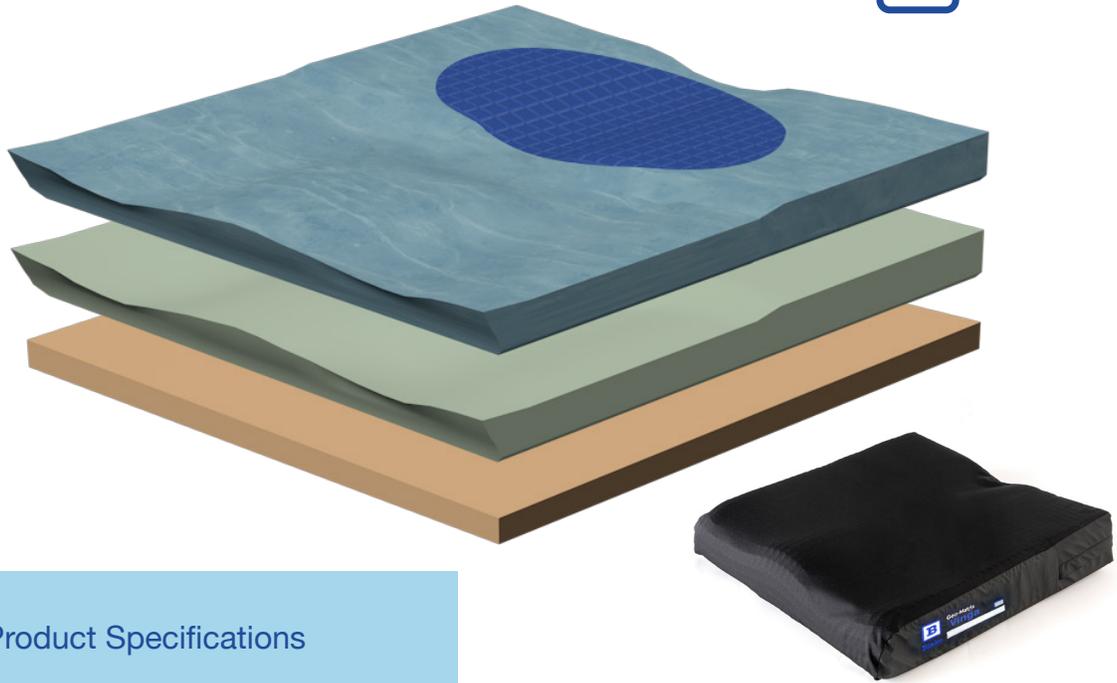

## Product Specifications

### CONSTRUCTION

Closed cell foam base and lateral stabilizing wedges. High resiliency foam for the base, and ischial shelf. Core is comprised of a contoured injected foam. The top layer features a gel infused memory foam at the front of the cushion and a zoned Pillow Top Geo-Matrix™ gel insert along the back half of the cushion.

### COVERS

Two outer cover options: Breathable Spacetex fabric with incontinent protection OR Elastatex soft touch outer incontinent cover.

Inner cover: Elastatex breathable, wipe-able, antibacterial, fluid proof cover.

### STANDARD SIZES

40x40, 40x45, 40x50,  
45x40, 45x45, 45x50,  
50x40, 50x45, 50x50,  
55x40, 55x45, 55x50,  
55x55 cm

### CUSTOM SIZES

All dimensions available

### PRODUCT CODE

GV(+size) eg. GV4545

### WEIGHT CAPACITY

Recommended Maximum 125kg

The Geo-Matrix™ Vinga cushion perfectly combines the pressure and shear reduction of Geo-Matrix™ gel with the plush comfort of gel infused memory foam.

The Vinga was designed to maximize comfort and protection while maintaining a stable posture. With its multi foam core and the firm ischial shelf the Vinga finds the balance between comfort and stability.

The Vinga is designed to be light-weight and maximize comfort, while maintaining stable posture.

Each Vinga cushion comes complete with either a breathable or soft touch incontinent cover, as well as a 4 way stretch protective inner cover.

Blake Medical esteemed distributor:



FOAM CUSHION  
**Geo-Matrix Vinga**




## WASHING INSTRUCTIONS

Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel and foam: Wipe down gel using any cleaner that does not contain bleach and allow to air dry.

**For more specific cleaning instructions please visit our website or user's manual.**

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam. All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

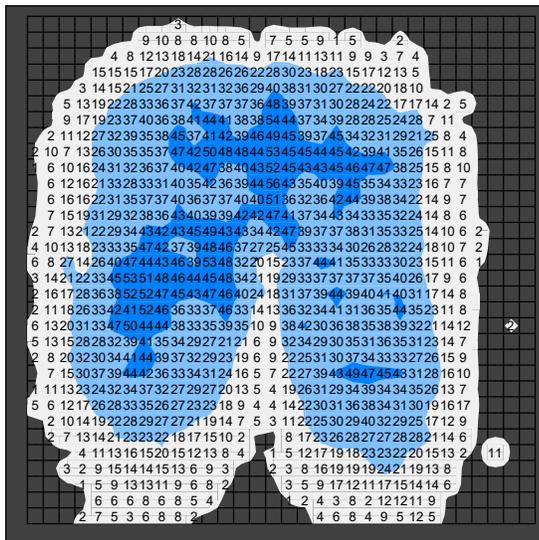

Client: Female 53kg

## ADDITIONAL BENEFITS

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management.
- Minimizes shear through column segmentation.
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold.
- Geo-Matrix™ gel is recyclable through Blake Medical.

## Advantages of Geo-Matrix™ Gel

 DURABLE	 COMFORTABLE	 DISINFECT-ABLE	 FRICTION & SHEAR REDUCING	 INCREASES SIT TIMES
 MICROCLIMATE MANAGEMENT	 PRESSURE REDISTRIBUTING	 VIBRATION DAMPENING	 WON'T PUNCTURE	

